



# SERVICE BULLETIN

---

## Health Imaging Products

---

Eastman Kodak Company, Health Imaging, Rochester, NY 14650

**SERVICE BULLETIN NO. 294**

**October 2000**

### ***3752 Kodak Min-R MAMMOGRAPHY PROCESSOR***

#### **Proper Installation of the JUMPER on TB3 in the Electrical Box**

The Kodak Min-R MAMMOGRAPHY PROCESSOR is shipped with a JUMPER on TB3.

For serial numbers 230001 through 230598, the processor has been shipped with the JUMPER in the 50 Hz position (TB3-7 to TB3-8).


For serial numbers 230458, 230460, 230461, 230463, 230467, 230468, 230471, 230472, 230482, 230483, 230485, 230489, 230491, 230494, 230496 through 230510 and 230599 and above, the processor is shipped with the JUMPER installed in the 60 Hz position (TB3-7 to TB3-9).

During installation of the processor, it is important to make sure the JUMPER is installed in the correct position, to match the power frequency at the site, or you will damage DRIVE MOTOR.

Also, all processors that have been previously installed, should be checked at the next service call or PM, to make sure the JUMPER is in the correct position:

- [1] Determine the frequency and voltage of the main power at the site and set the frequency and voltage of the PROCESSOR to match the power at the site:
- [2] If necessary, move the JUMPER to the correct position on TB3:

Frequency	Position
50 Hz	TB3-7 to TB3-8
60 Hz	TB3-7 to TB3-9

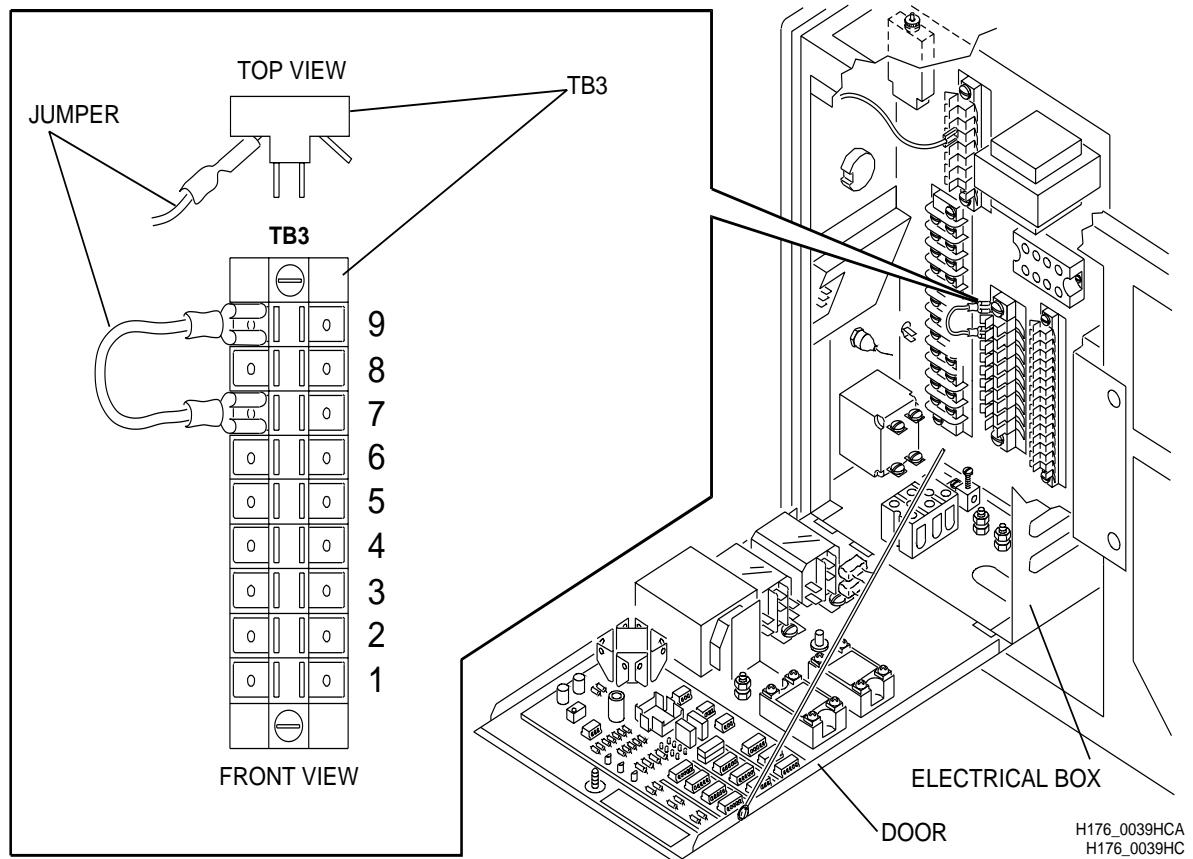

**CAUTION:** It is very important that there is only one (1) JUMPER installed on TB3. If necessary, you need to move the existing JUMPER to the correct location. **DO NOT** install an additional JUMPER or you will damage the DRIVE MOTOR.

**Debbie Trubatch, Service Engineer- Rochester**

**For More Information:**

**U.S. Distributors Contact:**

**Health Imaging  
Customer Support Operations  
1-800-336-4722**

**GCSS Personnel Contact:**

**Film Handling TAC  
U.S.: 1-800-822-1414  
Canada: 1-800-433-1414  
International: 716-724-4675**

*Kodak, X-Omat, and Min-R* are trademarks.

Pub No. 7E6503

SB\_#294.doc

Printed in USA

HEALTH IMAGING

EASTMAN KODAK COMPANY – ROCHESTER, NY 14650

